



EO-300

IPPB Device



Coaching function

Enabling the coaching function

The user can switch to Coaching screen to perform a therapy previously saved as preset by the medical personal.

NOTE: Coaching function is only available when a preset is selected.

From Home screen, click on  to display the coaching screen

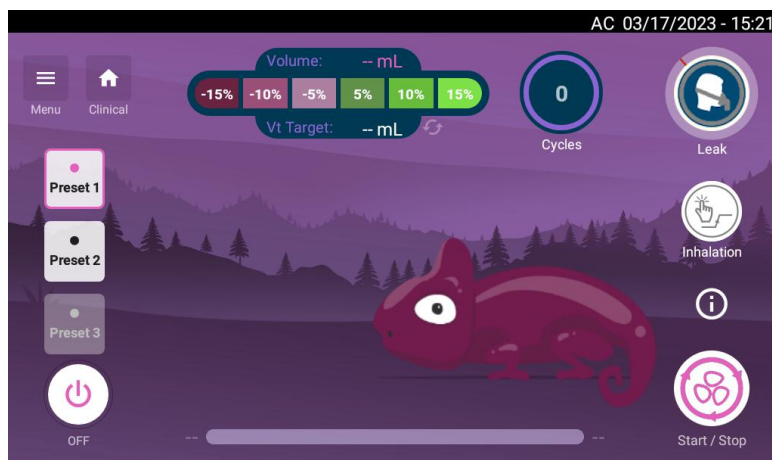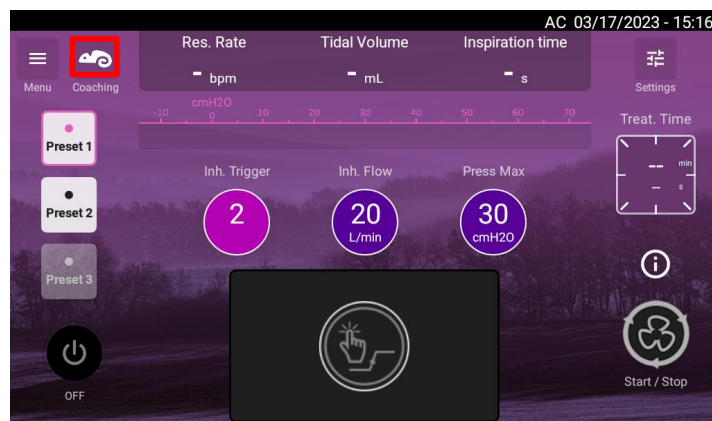

NOTE: The coaching function is Only available when a preset is selected


If the screen indicates no compatible preset is available, contact your service provider.

Changing preset mode from Coaching screen


1. Current preset /Activated preset

2. Available preset

3.Empty preset

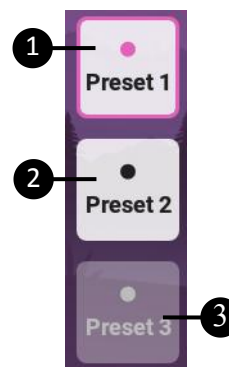

The IPPB Coaching Screen

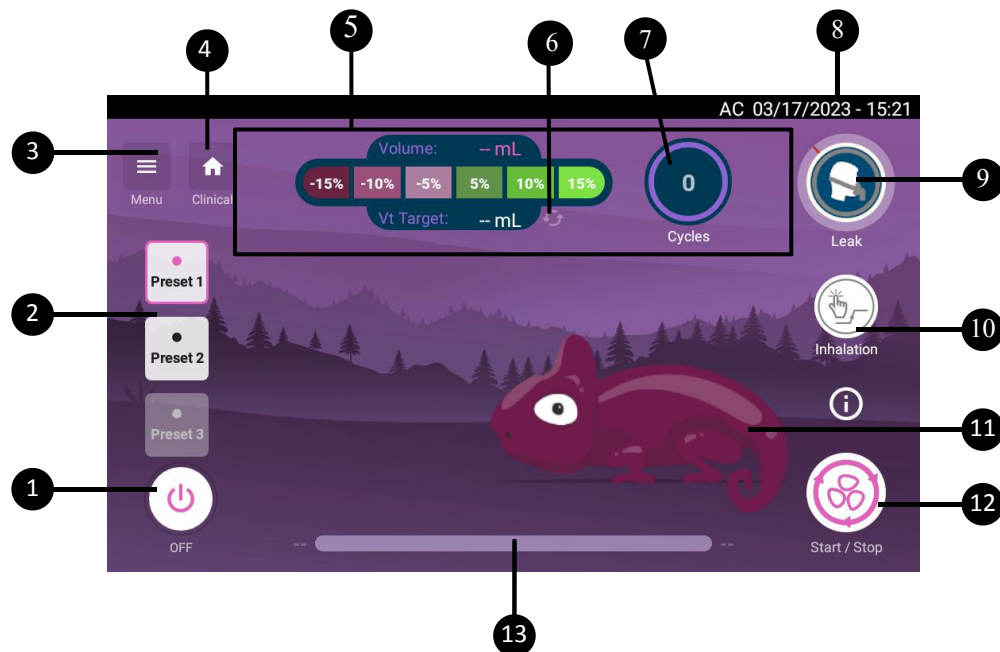

1. Power off button 2. Preset mode menu (1-3): Presets accessible to the patient when needed. 3. Menu button: Allows access to clinical home screen and clinical menus 4. Clinical screen button: Allow access to the clinical home screen 5. therapy monitoring: display the live value for the displayed monitoring during therapy. 6. Reset of the target volume previously defined 7. Number of breaths.	8. Date and Time with format YYYY/MM/DD (can be changed in Preferences menu) 9. Leak indicator 10. Trigger button: permit to manually trigger an inspiration during a therapy 11. Chameleon: Animation which will follow the therapy set up to help the patient 12. Start/Stop button: starts or stops therapy. 13. therapy Time Bar: display the elapsed and the remaining therapy time during therapy when a therapy time is set by the clinician
--	--

Performing a therapy with Coaching function

Select a preset and start the therapy. Following the therapy set up, inspiration can either be triggered by the patient or by using the dedicated button.

Volume shows the volume live evolution during the inspiration. This monitoring stays to the final inspiration volume until the next inspiration. The breath counter is incremented at the beginning of each inspiration. At the end of each inspiration, the patient interface leaks are evaluated and the

leak indicator  becomes green or orange. Orange will indicate that there is too much leak at the patient interface and that it should be better adjusted.

During inspiration, chameleon will inflates following the proportion of the target volume reached. It will progressively turn to green when reaching and exceeding the target volume.

For the 3 first inspirations:

- Chameleon does not move
- V Target stays undefined

After 3 inspirations:

- V target is defined. It is the mean of the volume of the 3 first inspirations. It stays unchanged until end of the therapy.

During exhalation, chameleon will gradually deflate until reaching the initial position and colour. It stays like this until next inspiration. The volume bargraph arrow will move following the proportion of the target volume reached. It will stay to the maximum level reached until next inspiration.

