



EO-70

Coaching function



User Guide

Enabling the coaching function

The user can switch to Coaching screen to perform a treatment previously saved as preset by the medical personal.


NOTE:	Coaching function is only available for preset with treatment in INEX-Auto mode or in IPPB mode.
--------------	--

From Home screen, click on  to display the Coaching function screen.

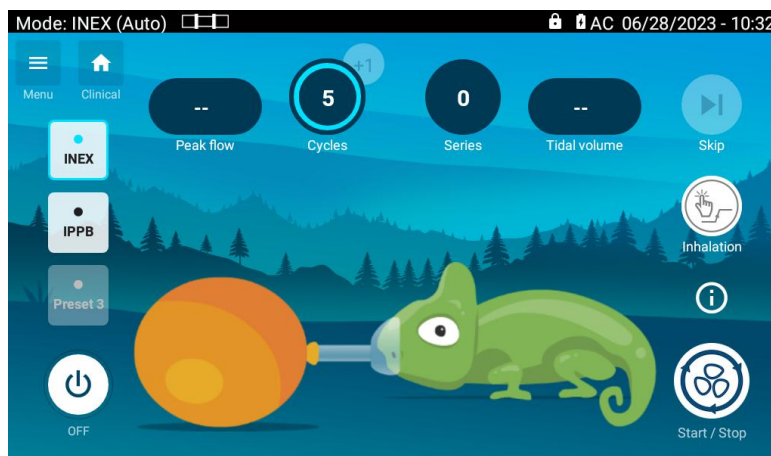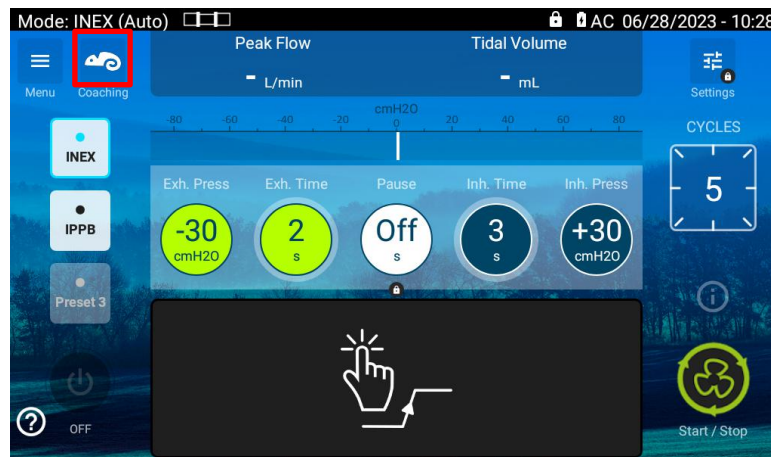

NOTE:	The coaching function is not available for all treatment modes. A message will indicate if something is missing to start the function.
	If the screen indicates no compatible preset is available, contact your service provider.

Changing preset mode from Coaching screen

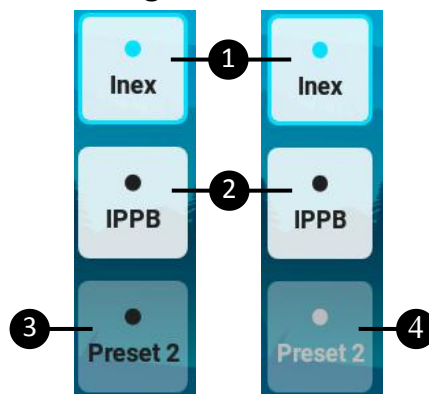


1. Current preset / Activated preset
2. Available preset
3. Inactive preset (not available for Coaching)
4. Empty preset

To change the preset, click on the preset you want to switch on.

If the treatment is in progress, it is not possible to change the preset.

The INEX Coaching Screen

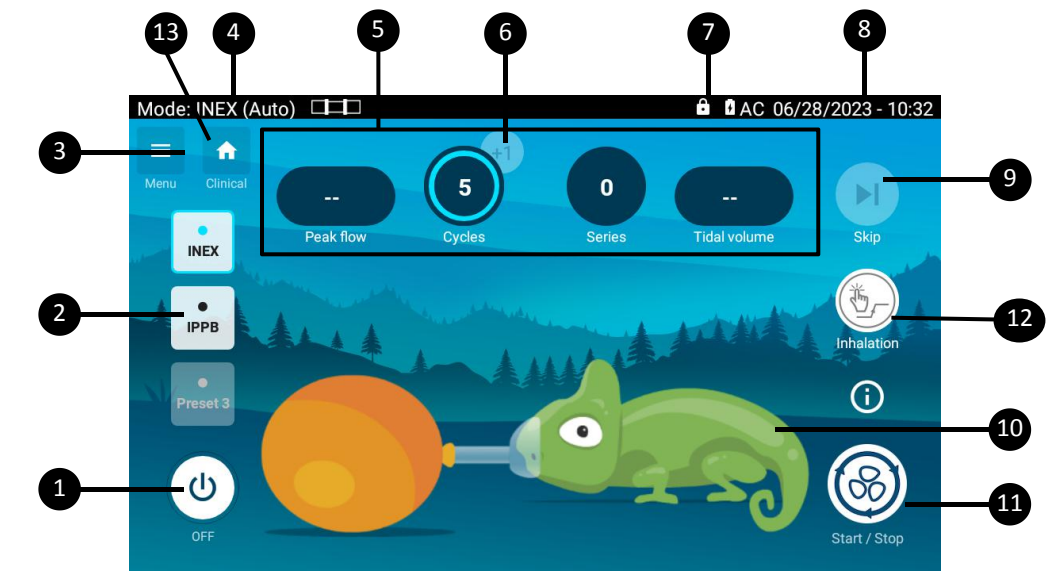


1. Power off button 2. Preset mode menu (1-3): Presets installed by the clinician and accessible to the patient when needed. 3. Menu button: Allows access to patient screen and clinical menus 4. Treatment mode indicator: INEX (manual/Auto) or IPPB 5. Treatment monitoring: display the live value for the displayed monitoring during treatment. 6. Cycle addition button: permit to add a cycle to the ongoing treatment if necessary	7. Battery life indicators: Indicates the level of charge left in the battery or whether the battery is charging. 8. Date: indicates the date on YYYY/MM/DD format. Time: Indicates the time on 24hr clock. This can be set up and changed from the Preferences menu. 9. Skip button: permit to stop the INEX serie on going. 10. Chameleon: Animation which will follow the treatment set up to help the patient 11. Start/Stop button: starts or stops treatment. 12. Trigger an inspiration during the treatment 13. Clinical screen shortcut
---	--


Performing an INEX treatment with Coaching function

Select a preset with an INEX Auto treatment and start the treatment. Chameleon and balloon animation will follow the treatment phases:

- During inspiration, chameleon inflates, and balloon deflates

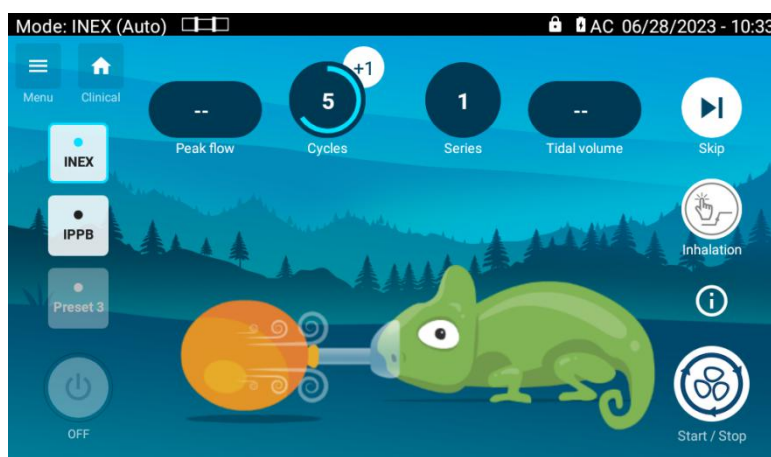

- During exhalation, chameleon deflates, and balloon inflates

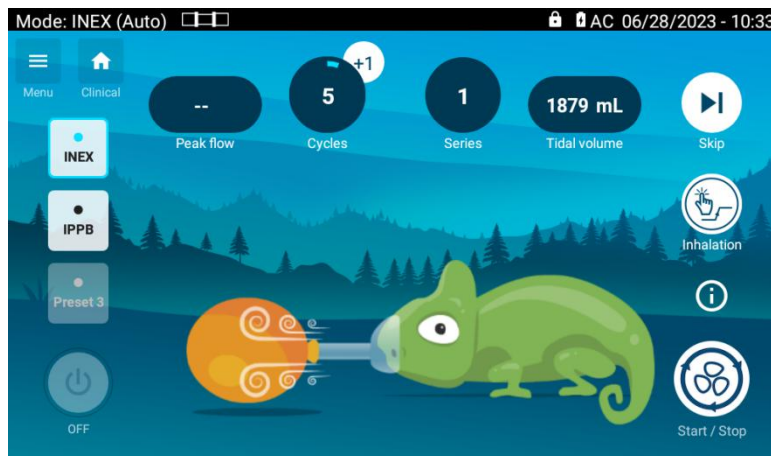


- During pause, chameleon and balloon do not move.

The monitorings:


- The peak flow is updated after each exhalation
- The cycle decoupler decreases after each exhalation. It can be increased with the «+1» button during the last cycle
- Series counter is incremented at each start of treatment. It is reset when changing preset or shutting down the devices
- Tidal volume is updated after each inspiration

At the end of the serie (or when stopping the treatment using the skip button), a pop up will appears to ask the patient if a new serie will be performed or not.

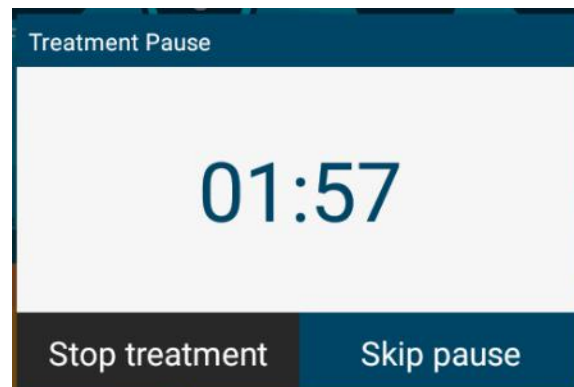

The chronometer is indicative. Click on «skip pause» to restart a serie or on «Stop treatment» to stop the treatment.

The IPPB Coaching Screen

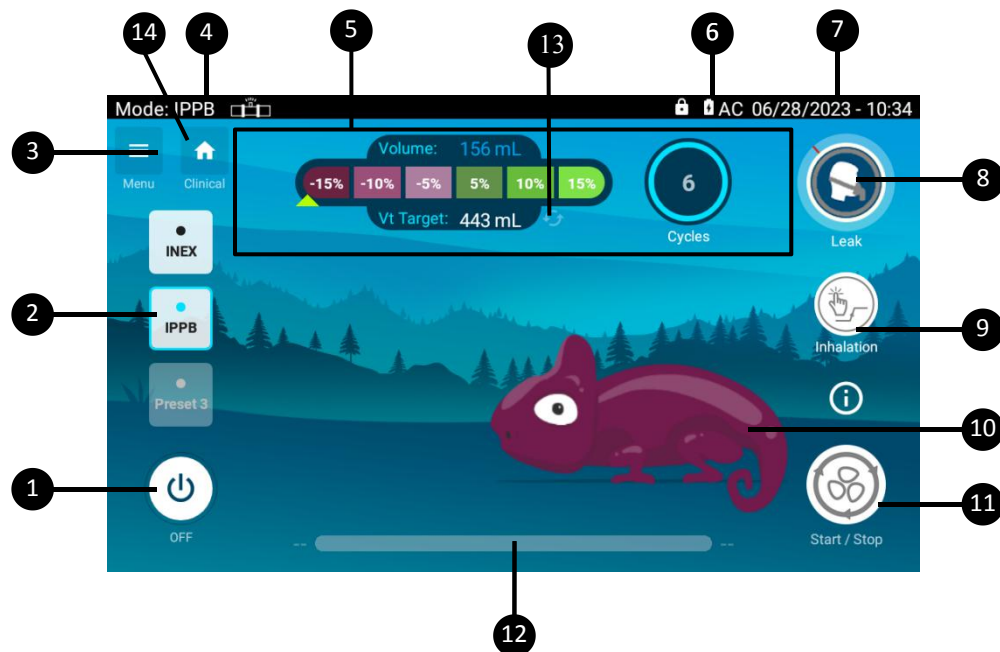

1. Power off button 2. Preset mode menu (1-3): Presets installed by the clinician and accessible to the patient when needed. 3. Menu button: Allows access to patient screen and clinical menus 4. Treatment mode indicator: INEX (manual/Auto) or IPPB 5. Treatment monitoring: display the live value for the displayed monitoring during treatment. 6. Number of breaths 7. Date and Time with format YYYY/MM/DD (can be changed in Preferences menu)	8. Leak indicator 9. Trigger button: permit to manually trigger an inspiration during a treatment 10. Chameleon: Animation which will follow the treatment set up to help the patient 11. Start/Stop button: starts or stops treatment. 12. Treatment Time Bar: display the elapsed and the remaining treatment time during treatment when a treatment time is set by the clinician 13. Reset of the target volume previously defined 14. Clinical screen shortcut
--	--

Performing an IPPB treatment with Coaching function

Select a preset with IPPB treatment and start the treatment. Following the treatment set up, inspiration can either be triggered by the patient or by using the dedicated button.

Volume shows the volume live evolution during the inspiration. This monitoring stays to the final inspiration volume until the next inspiration. The breath counter is incremented at the beginning of each inspiration. At the end of each inspiration, the patient interface leaks are evaluated and the

leak indicator  becomes green or orange. Orange will indicate that there is too much leak at the patient interface and that it should be better adjusted.

During inspiration, chameleon will inflate following the proportion of the target volume reached. It will progressively turn to green when reaching and exceeding the target volume.

For the 3 first inspirations:

- Chameleon does not move
- V Target stays undefined

After 3 inspirations:

- V target is defined. It is the mean of the volume of the 3 first inspirations. It stays unchanged until end of the treatment.

During exhalation, chameleon will gradually deflate until reaching the initial position and colour. It stays like this until next inspiration. The volume bargraph arrow will move following the proportion of the target volume reached. It will stay to the maximum level reached until next inspiration.

